Información general:

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All for Healthcare

Cara Trasera ampliada

EN • RS6 Wrist blood pressure monitor. Clinically validated • Cuff wrapping guide • Irregular Heartbeat Detector • For wrist circumferences 13,5-21,5 cm • Hardcase • Batteries Included

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Manual de Instrucciones

Introduction

Thank you for purchasing the OMRON® Wrist Blood Pressure Monitor.

The Wrist Blood Pressure Monitor is a compact and easy-to-use blood pressure monitor, operating on the oscillometric principle. It measures your blood pressure and pulse rate simply and quickly. For comfortable and accurate readings, the use of a pressure cuff with or without the device uses an advanced "intelligent" technology.

Important Safety Information

Consult your doctor prior to using a pressure monitor, and if diagnosed with any condition of the heart or blood vessels. Please read the entire manual before using the unit.

Warnings:

- This device is not intended for use by pregnant women.
- Do not use the device if you are pregnant or if you are taking any medication that affects blood pressure.
- Do not use the device if you are taking any medication that affects blood pressure.

1. Overview

Main unit:

A. Main unit
B. Battery compartment
C. Display
D. Inflation/Deflation buttons
E. Power button
F. On/Off button
G. Memory button
H. Cuff mapping button

Display:

1. Memory symbol
2. System blood pressure
3. Diastolic blood pressure
4. Memory symbol
5. Inflation/Deflation buttons
6. Power button
7. On/Off button
8. Memory button
9. Cuff mapping button
10. Cuff mapping button
11. Cuff mapping button
12. Cuff mapping button
13. Cuff mapping button
14. Cuff mapping button
15. Cuff mapping button
16. Cuff mapping button
17. Cuff mapping button
18. Cuff mapping button
19. Cuff mapping button
20. Cuff mapping button

2. Preparation

1. Remove the battery cover.


2. Insert two 1.5V "AAA" alkaline (LR03) batteries as indicated in the battery compartment and then replace the battery cover.

3. Setting the Date and Time

1. Press the **ON/OFF** button.

2. Set the monitor to the correct date and time before taking a measurement for the first time.

3. Press the **START/STOP** button to store the setting.

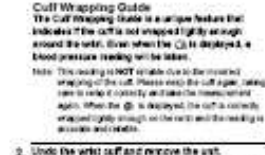3. Using the Unit

3.1 Applying the Wrist Cuff

1. Put your arm through the cuff loop.

2. Wrap the wrist cuff around your wrist.

3.2 How to Sit Correctly

1. Take a measurement, you need to be relaxed and comfortable.

2. Take a measurement, you need to be relaxed and comfortable.

3.3 Taking a Reading

1. Press the **START/STOP** button.

2. Press the **START/STOP** button to store the setting.

Important Safety Information

Consult your doctor prior to using a pressure monitor, and if diagnosed with any condition of the heart or blood vessels. Please read the entire manual before using the unit.

Warnings:

- This device is not intended for use by pregnant women.
- Do not use the device if you are pregnant or if you are taking any medication that affects blood pressure.
- Do not use the device if you are taking any medication that affects blood pressure.

General Precautions:

- Do not use the device if you are pregnant or if you are taking any medication that affects blood pressure.
- Do not use the device if you are taking any medication that affects blood pressure.
- Do not use the device if you are taking any medication that affects blood pressure.

2. Preparation

1. Remove the battery cover.

2. Insert two 1.5V "AAA" alkaline (LR03) batteries as indicated in the battery compartment and then replace the battery cover.

3. Setting the Date and Time

1. Press the **ON/OFF** button.

2. Set the monitor to the correct date and time before taking a measurement for the first time.

3. Press the **START/STOP** button to store the setting.

3. Using the Unit

3.1 Applying the Wrist Cuff

1. Put your arm through the cuff loop.

2. Wrap the wrist cuff around your wrist.

3.2 How to Sit Correctly

1. Take a measurement, you need to be relaxed and comfortable.

2. Take a measurement, you need to be relaxed and comfortable.

3.3 Taking a Reading

1. Press the **START/STOP** button.

2. Press the **START/STOP** button to store the setting.

Cuff Wrapping Guide

The cuff wrapping guide is a unique feature that indicates if the cuff is not wrapped tightly enough around the wrist. Once when the **ON/OFF** displays a blood pressure reading will be taken.

Important:

- Do not use the device if you are pregnant or if you are taking any medication that affects blood pressure.
- Do not use the device if you are taking any medication that affects blood pressure.
- Do not use the device if you are taking any medication that affects blood pressure.